




SAY SKIN



SAY SKIN

saying beauty

CONTENTS

- COMPANY PROFILE
- OVERVIEW
- DESCRIPTION OF PRODUCTS
- REFERENCE
- CORE TECHNOLOGIES
- HOW IT WORKS
- COMPARISON TABLE
- CLINICAL TRIAL
- HOW TO USE
- PATENTS & CERTIFICATES

Company Profile

SAY SKIN



Company Name	WITHNIX Co., Ltd.
Date of Establishment	October 1st, 2003
Head Office	593-14, Opo-ro, Opo-eup, Gwangju-si, Gyeonggi-do, Republic of Korea
Subsidiary	WITHNIX Shenzhen Co., Ltd. (China)
Home Page	www.withnix.co.kr en.say-skin.com



WITHNIX CO.,LTD is a subsidiary of major pharmaceutical company in Korea, Dongkook pharmaceutical.

We shares the core management philosophy of its parent company,
which is dedicated to ensuring the healthy and happy lives of its customers.

AURORA⁺



AURORA-*lite*



Description

AURORA⁺



BRAND	SAY SKIN
CATEGORY / PRODUCT	Beauty Device / AURORA+
BENEFITS FOR SKIN	<p>* PROMOTION ABSORPTION OF COSMETICS Boosts the absorption of effective ingredients of skin are products to provide strong protection for the dermis.</p> <p>* TROUBLE CARE Solves problems of skin trouble and irritation caused by Exposure to pollutants, such as fine dust, once and for all.</p> <p>* WRINKLE TREATMENT & SKIN TIGHTENING Helps tighten skin and reduces wrinkles.</p> <p>* WHITENING AND SKIN TONE UP Improves dull complexion and skin blemishes and provides A strong skin whitening effect.</p>
TECHNOLOGIES	Plasma, Microcurrent
MATERIALS	[Body] PC+ABS / [Heads] Zirconia ceramic
SIZE / WEIGHT	47.0(W) * 184.5(L) * 56.1(H)mm / 185g
POWER TYPE	DC 5V, 1A (USB-C type)
BATTERY	Li-Polymer Battery 1,000mAh
CERTIFICATION	KC, CE, FCC
COUNTRY OF MANUFACTURE	Republic of KOREA
COLOR	WHITE, BLACK

Description

AURORA-*lite*



BRAND	SAY SKIN
CATEGORY / PRODUCT	Beauty Device / AURORA-lite
BENEFITS FOR SKIN	<p>* PROMOTION ABSORPTION OF COSMETICS Boosts the absorption of effective ingredients of skin are products to provide strong protection for the dermis.</p> <p>* TROUBLE CARE Solves problems of skin trouble and irritation caused by Exposure to pollutants, such as fine dust, once and for all.</p> <p>* WRINKLE TREATMENT & SKIN TIGHTENING Helps tighten skin and reduces wrinkles.</p> <p>* WHITENING AND SKIN TONE UP Improves dull complexion and skin blemishes and provides a strong skin whitening effect.</p>
TECHNOLOGIES	Plasma
MATERIALS	[Body] PC+ABS, Aluminum / [Heads] Zirconia ceramic
SIZE / WEIGHT	30.3(W) * 155(L) * 52(H)mm / 103g
POWER TYPE	DC 5V, 1A (USB-C type)
BATTERY	Li-Polymer Battery 700mAh
CERTIFICATION	KC, CE, FCC, CVC
COUNTRY OF MANUFACTURE	Republic of KOREA
COLOR	WHITE, PINK

Aesthetic Test for improving Skin trouble


*Using AURORA+ with the serum for 10~13minutes everyday



- **Target** : A boy in his 10's / Forehead and Cheek
- **Solution** : 7 times with the solution
- **Result** : Comparison of changes before and after caring with the solution on the face

Aesthetic Test for improving Skin trouble

*Using AURORA+ with the serum for 10~13minutes everyday

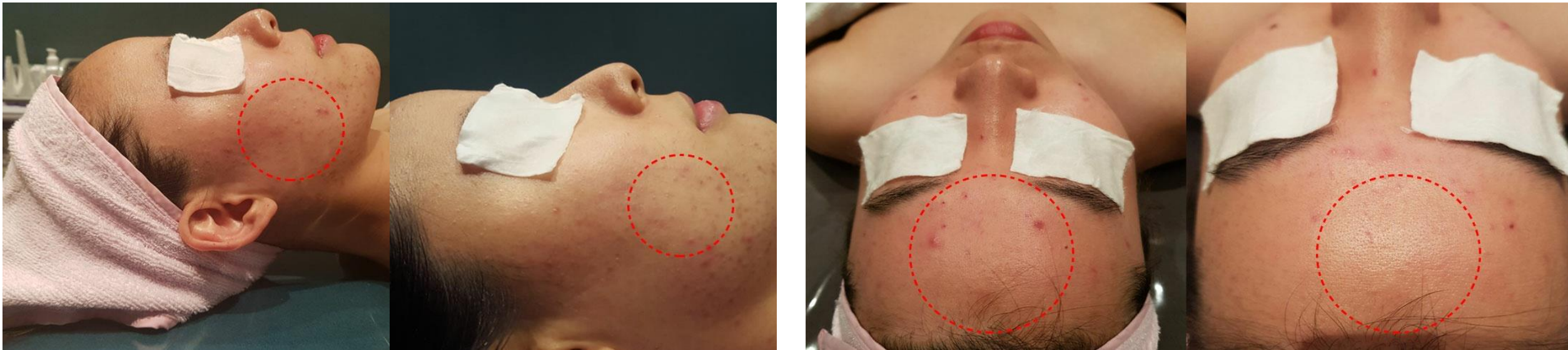

- **Target** : A woman in her 20's / Forehead and Cheek
- **Solution** : 12 times with the solution
- **Result** : Comparison of changes before and after caring with the solution on the face

Aesthetic Test for improving Skin trouble

*Using AURORA+ with the serum for 10~13minutes everyday



- **Target** : A woman in her early 30's / Forehead and Cheek
- **Solution** : 13 times with the solution
- **Result** : Comparison of changes before and after caring with the solution on the face

Self-Test for absorption ability of cosmetics

60 seconds after dropped a drop of serum into the back of right hand, **without** cared AURORA-lite.



VS

3 seconds after dropped a drop of serum into the back of left hand, **cared with** AURORA-lite for 2 minutes.



* **Result** : The serum is immediately absorbed into the skin cared with AURORA-lite.

PLASMA



**Cutting-Edge, the 4th material,
Proven harmless to the human body.**

- Promotion of skin regeneration.
- Strong sterilization effect.
- Promotion of cosmetic absorption.

MICROCURRENT



A Mega-trend of beautifying device.

- Promotion of Skin elasticity and pore reduction.
- Promotion of physical metabolism.


ZIRCONIA




A Safe and Eco-friendly Ceramic.

- A high-strengthened, precise and eco-friendly ceramic to emphasize safety.

How it works

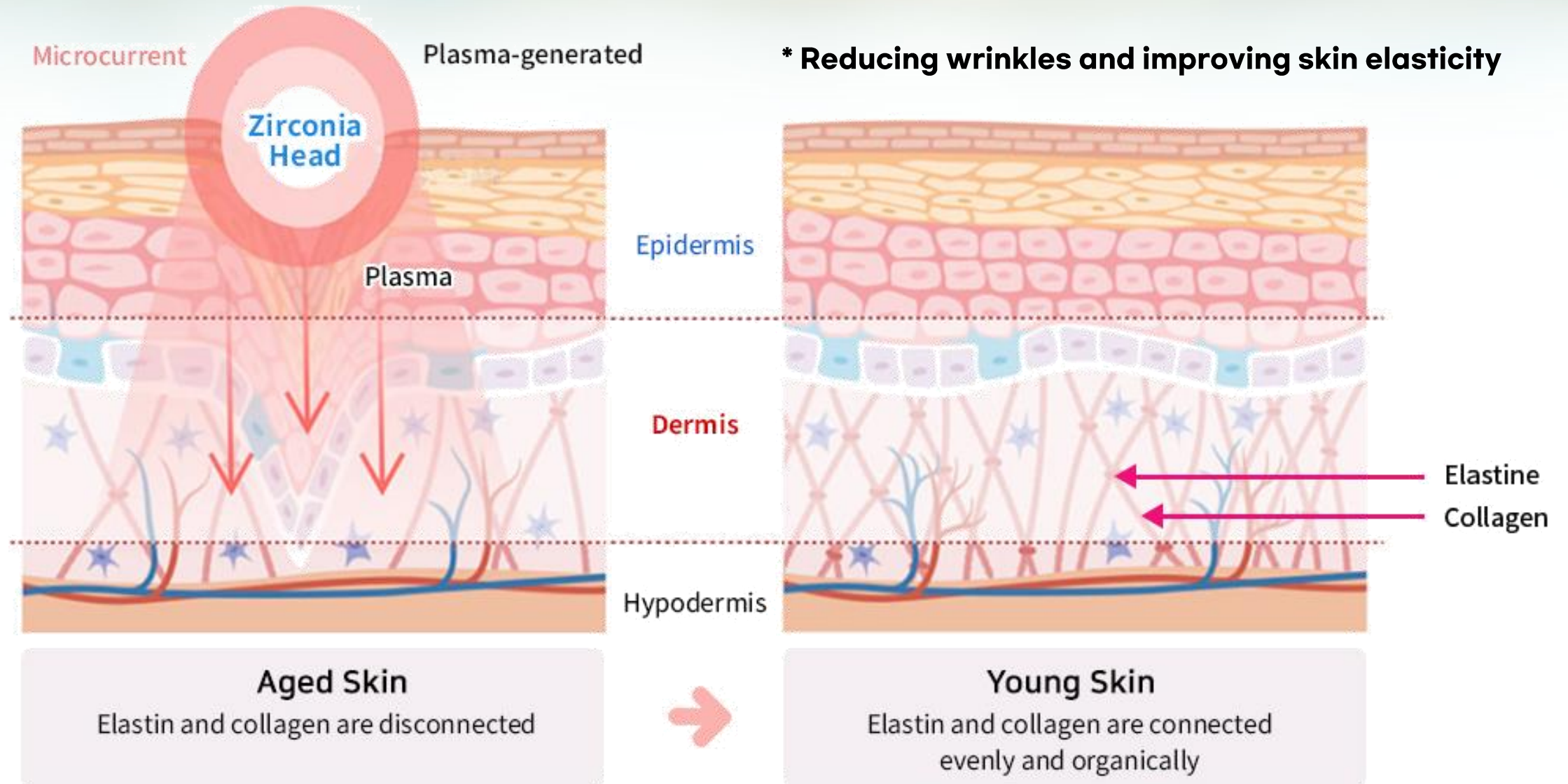

How it works

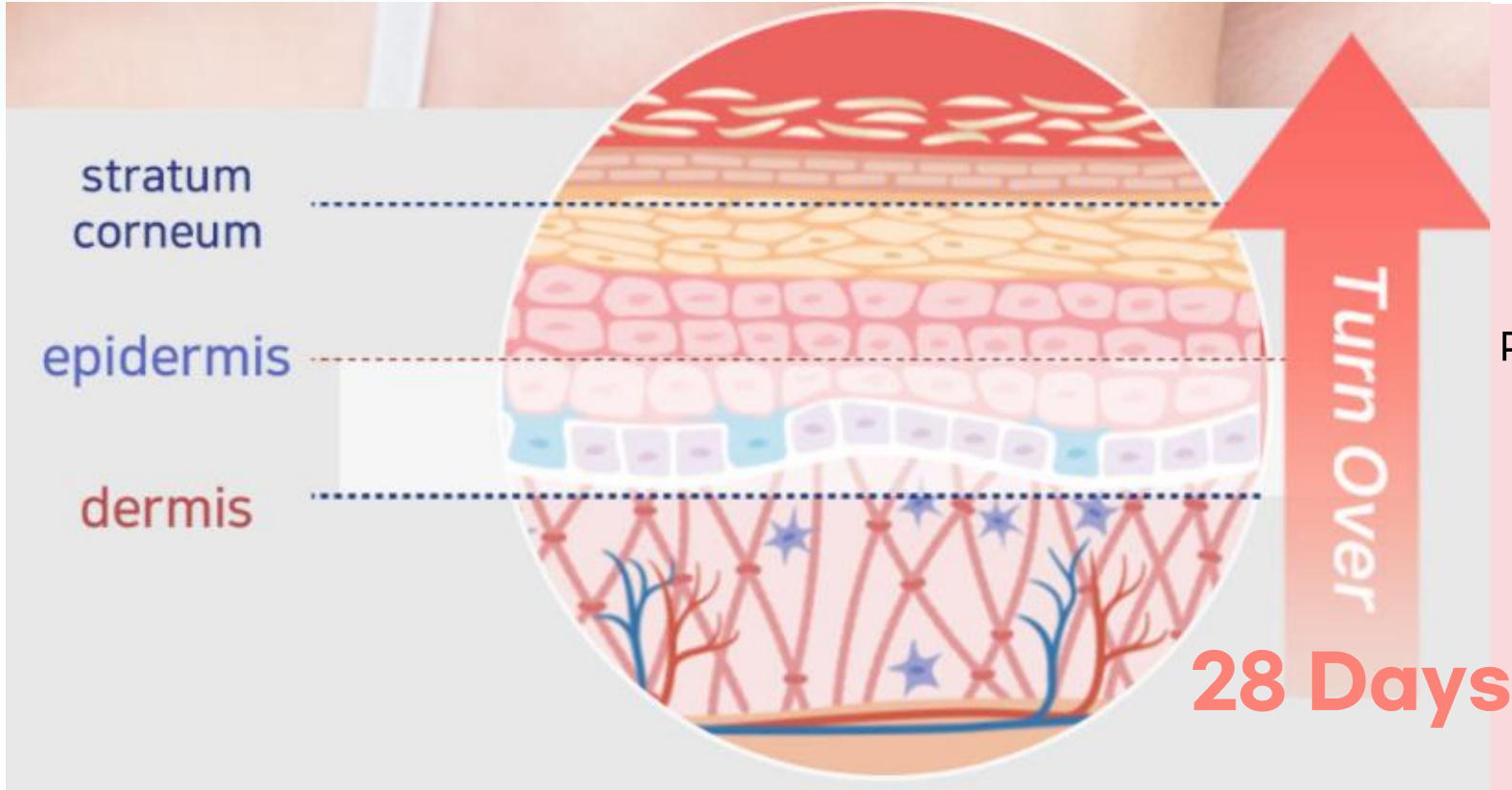

*** Whitening and improving skin texture**


Plasma causes the desquamation which eliminates keratin mass, dead keratin cells and various kinds of contaminant.

How it works

Before



After

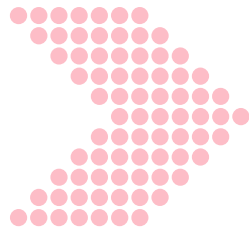

Source : Non-Thermal Plasma for Acne and Aesthetic Skin Improvement


*** Making skin troubles get better**

A hydroxyl group (OH radical), an oxygen-anionic natural substance produced when plasma is generated, has the strongest sterilizing power on the earth.

Sterilizes, deodorizes and decomposed by directly engaging in contaminants in air and water.

Sterilizing rate 2,000 times and 180 times faster than ozone and solar UV ray

Cell Membrane

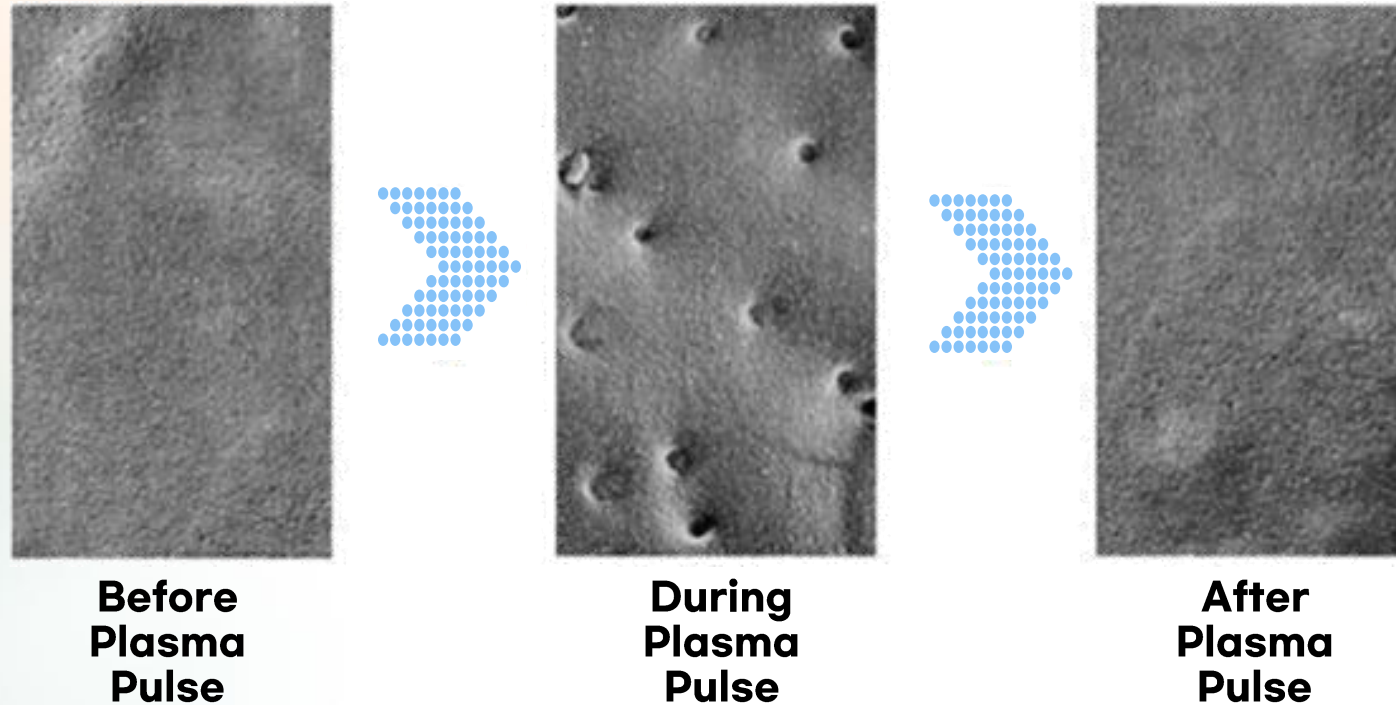







* Improving the skin's ability to absorb cosmetic active ingredients

“Plasmaporation = Plasma + Electroporation”

**Electroporation: the creation of transient pores in the cell membrane, through the application of a high-voltage electrical pulse.*

Comparison with Competitors

Model Name	SAY SKIN AURORA ⁺		SAY SKIN AURORA- <i>lite</i>	NuFace Trinity	MEDICUBE Age-r Booster Pro	LG Pra.L LED Mask	Dual Sonic Maximum
Pictures							
Country	WITHNIX / KOREA			USA	KOREA	KOREA	KOREA
Technology	PLASMA, MICROCURRENT, PLASMA-PORATION, GUASHA		PLASMA, PLASMA-PORATION	MICROCURRENT, GUASHA	MICROCURRENT, ELECTROPORATION, LED	LED	HIFU, RF, MICROCURRENT, LED
Effectiveness	Promotion Cosmetics Absorption, Elasticity, Whitening, Trouble Care (e.g., atopic dermatitis/acne) etc			Effectiveness on blood circulation, improvement of elasticity	Whitening, Promotion cosmetics absorption, Collagen activeness	Anti-Aging	Whitening, Promotion cosmetics absorption, Improvement of Elasticity
MSRP	\$ 250		\$ 145	\$ 260	\$ 290	\$ 510	\$ 1,100

Comparison with Competitors

(Plasma Devices)

Source : Withnix research and development department

1.

Division	Pla Beau	AURORA ⁺
Number of pulse per second	3,609	4,200
Number of pulse per second / Rate	85%	100%
Vrms of each pulse	545	590
Vrms of each pulse / Rate	92%	100%
Total PLASMA generation energy ratio	79%	100%
	100%	126%

* While the quantity or intensity of the plasma generated cannot be directly measured, the energy required for plasma generation can be considered to vary according to the potential applied to the head surface. Therefore, it is deemed valid to evaluate by comparing the electrical energy values measured on the head surface.

CONCLUSION

AURORA+ supplies 1.259 times stronger plasma generation energy compared to Pla Beau.

2.

Division	Plasma Shower	AURORA ⁺
Number of pulse per second	5,152.5	4,200
Number of pulse per second / Rate	122%	100%
Vrms of each pulse	350	590
Vrms of each pulse / Rate	59%	100%
Total PLASMA generation energy ratio	72%	100%
	100%	137%

CONCLUSION

AURORA+ supplies 1.373 times stronger plasma generation energy compared to Plasma Shower.



SAY SKIN

Clinical Trial

AURORA+ TEST

- Clinical institution : Korea Institute of Dermatological Sciences
- Clinical period : 6 weeks
- Test subject : 21 people, women aged 40~69 years old
- Evaluation item : Lifting of eye bag area under eyes, nasolabial area lifting, eye lifting, skin elasticity, skin density, skin tone



AURORA LITE TEST

- Clinical Institution : H&Bio SRC
- Clinical Period : June 7th ~ 21st, 2021 (2 weeks)
- Test subject : Dermalab Full Skin 3cm x 3cm / Acetyl-Hexapeptide-8(Lys-FITC)
- Evaluation Item : The ability of absorption cosmetics active ingredients
- Clinical institution : OATC
- Clinical period : 2 weeks
- Clinical subject : 22 people, women aged 30~60 years old
- Evaluation item : The effect of skin elasticity for AURORA-lite with the ampoule


The Effect of Skin Lifting

REF. Color of improved skin : Yellow (the best) > Green > Blue > Purple (the worst)

Under Eye Bags

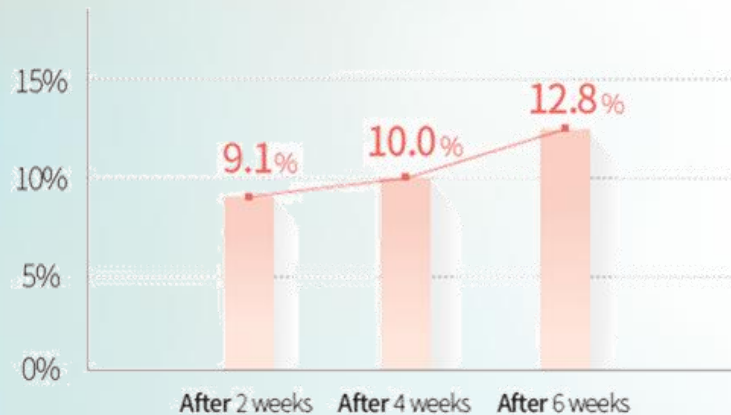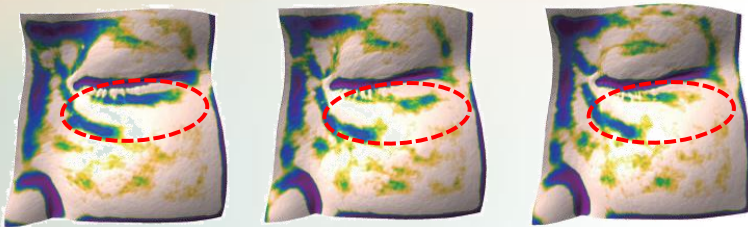


Nasolabial Folds

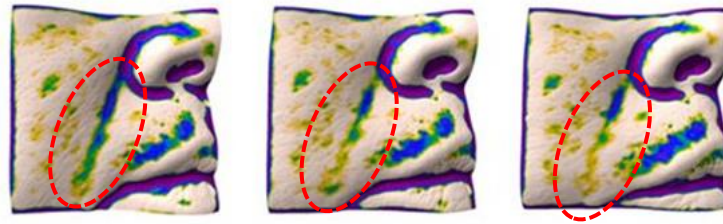


Before use Immediately after use 6 weeks after use



Crow's Feet

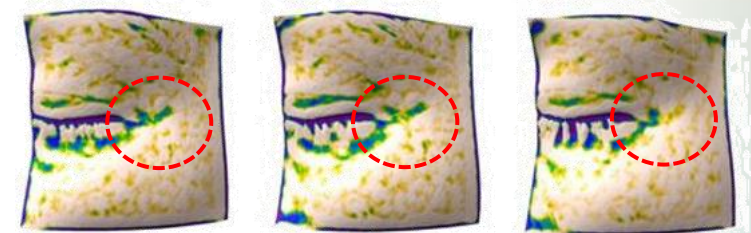


Before use Immediately after use 6 weeks after use

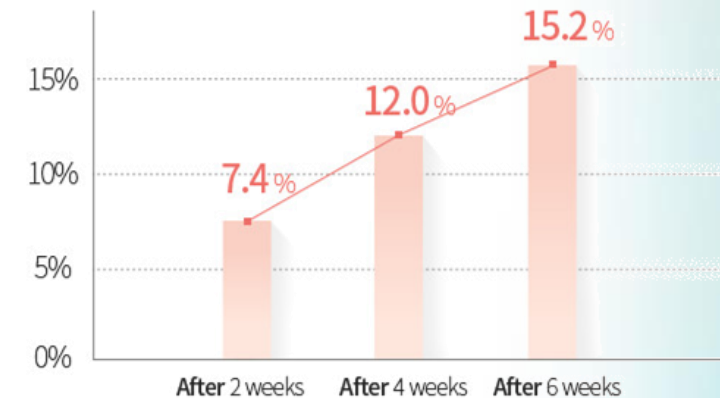


The Effect of

Skin Density

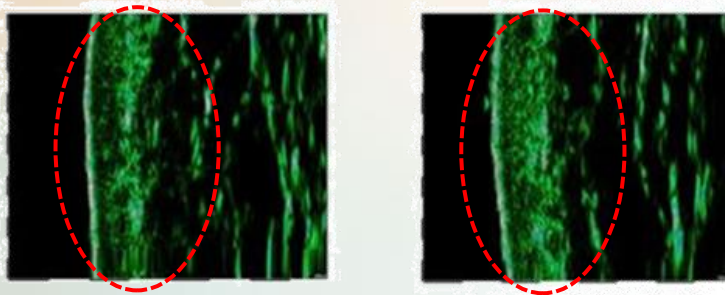


Before use


6 weeks after use

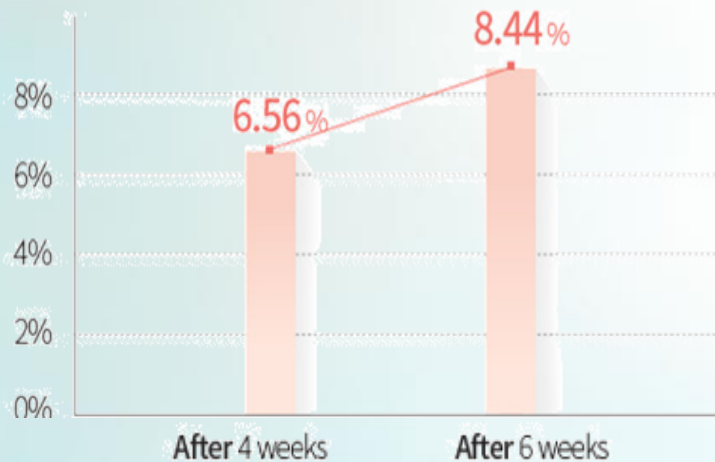

Skin Elasticity


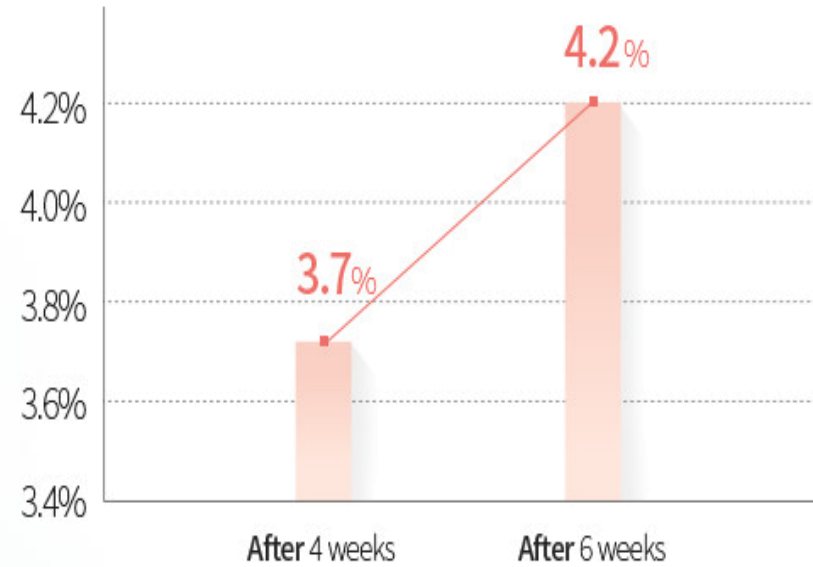


Skin Tone

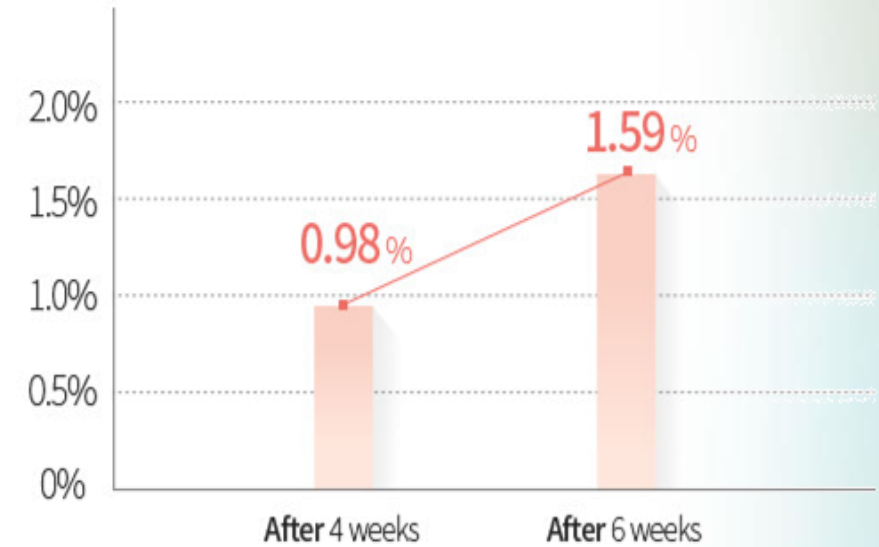

The ability of absorption of cosmetics active ingredients

* Comparison of skin absorption ex vivo human skin tissue

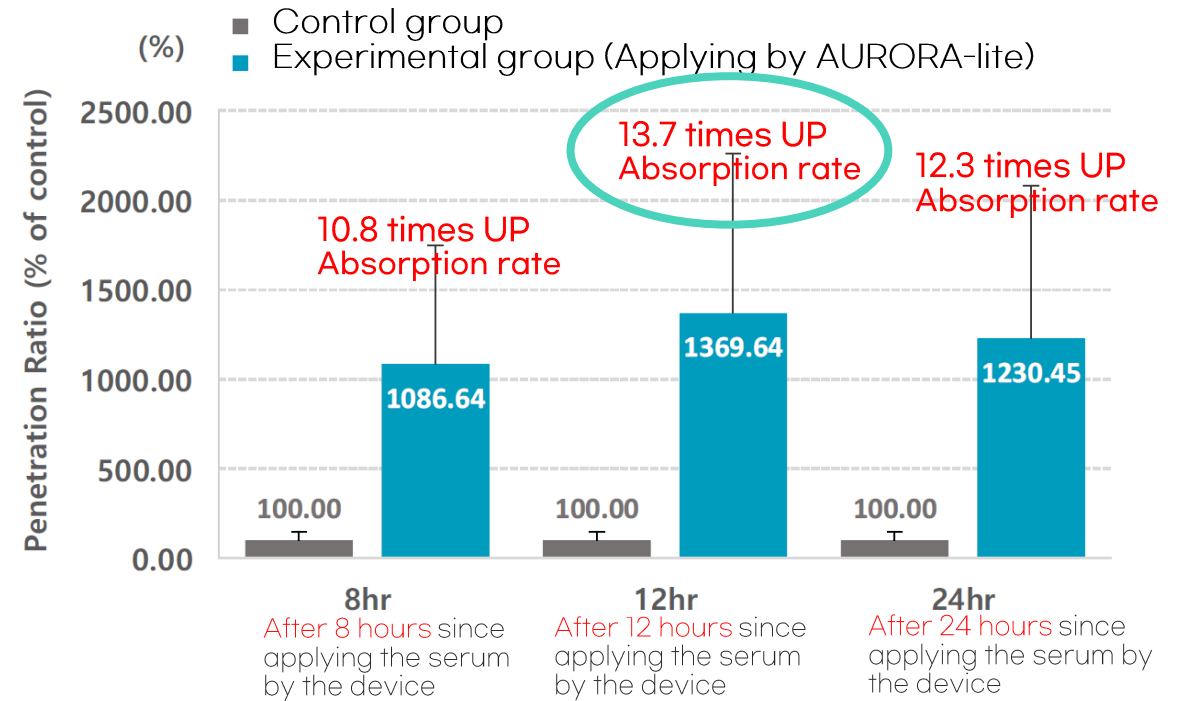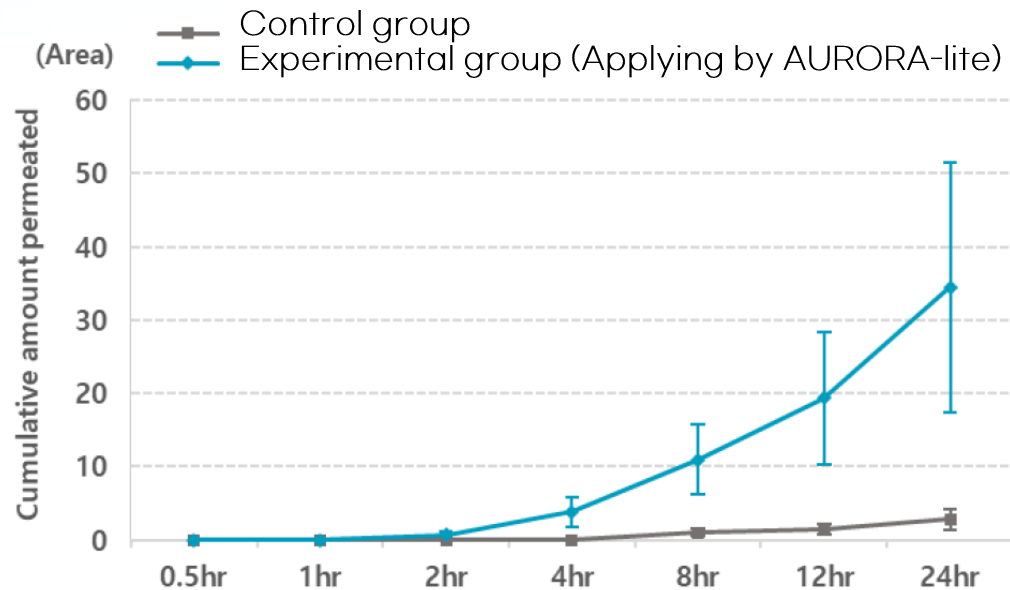


* More being absorbed the active ingredients of the serum into skin by the device than applying the serum by hand as below.


The ability of absorption of cosmetics active ingredients

After 1 hour

The group by hand

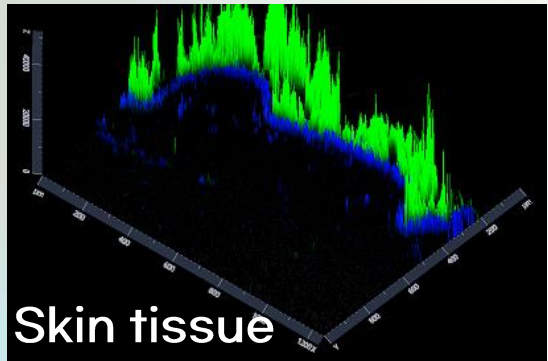

The group by AURORA-lite

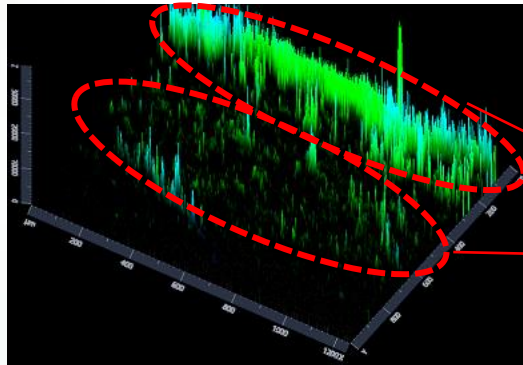


VS

Epidermis


Dermal Layer

After 24 hours

The group by hand

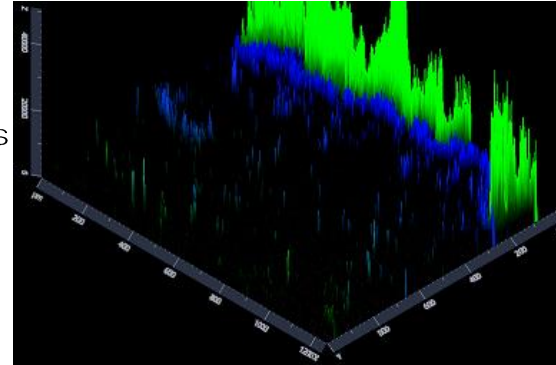

The group by AURORA-lite

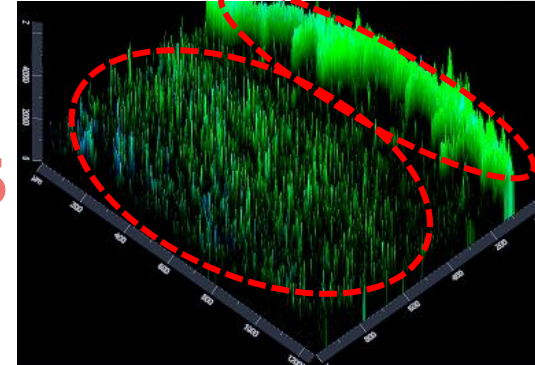


VS

REF. The substance in green is like a cosmetic active ingredient.

The Effect of Skin Elasticity

* Comparison of skin elasticity improvement

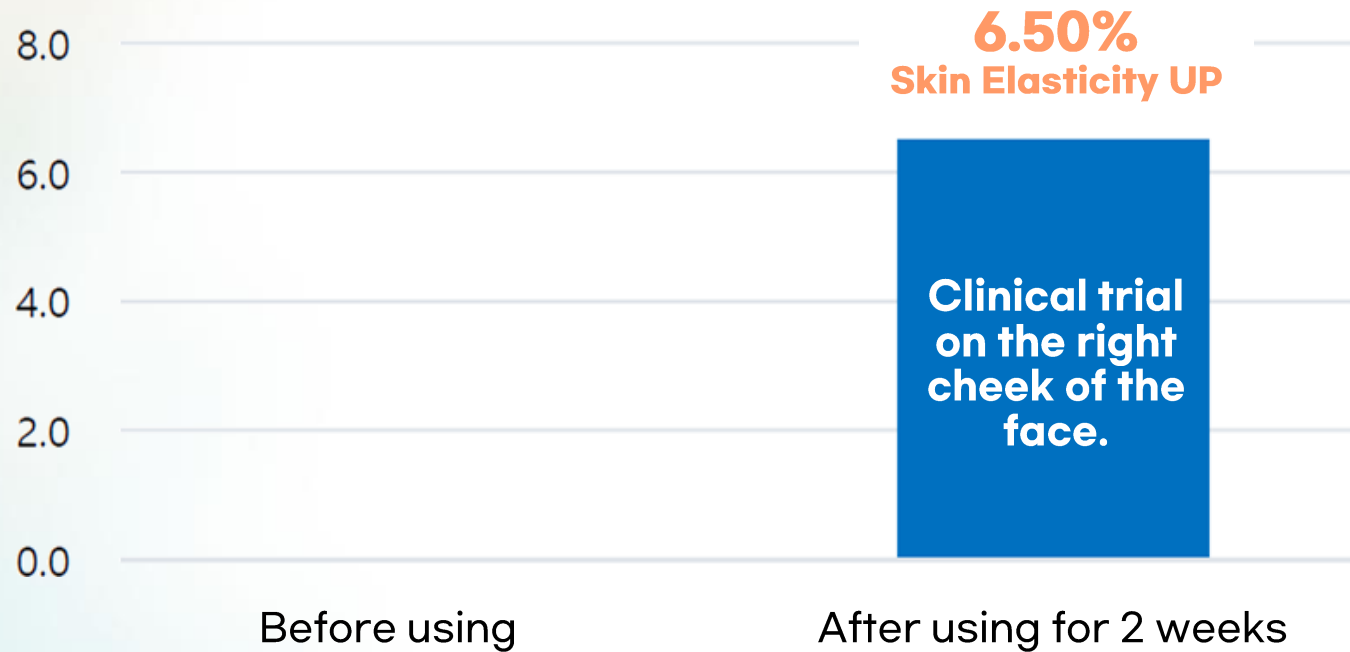

Post-clinical survey

SAY SKIN



* The results of the post-clinical survey showed the **generally high satisfaction** level with this product.

Material & Design

SAY SKIN

HEAD

- ✓ Using **Zirconia**, a harmless and eco-friendly material.
- ✓ **A cool feeling** reduces warmth during skin massage.
- ✓ Optimized structure of spacing and angle **between heads and head ball size** to care for curved areas such as **jawline, bridge of the nose, cheekbones.**

BODY

- ✓ Safe **PC+ABS** material.
- ✓ **UV coding** treatment prevents surface deterioration by cosmetics.
- ✓ Optimized design of long body with **a streamlined structure and grip to care for all parts of the body** such as nape, the back of the neck and the clavicle.

Operation Button

- ✓ With the buttons and displays located on the front, convenient UI design to operate with one thumb while holding the product.



Area-Specific Intensive care

- **Eye Area Care**

Close your eyes and use **one or two headballs**, moving them **back and forth** under your eyes to the eye bags.

- **Neck Care**

Avoiding the **thyroid area** in the front of the neck, please use the headballs in an upward motion from **bottom to top**.

- **Forehead Care**


Please use it in a repeating motion from the **eyebrow line** to the **hairline**, moving upward from bottom to top.

- **Cheek Care**

Please care for the headballs by moving from **left to right** and from **inside to outside** in a repetitive motion.

- **Jawline Care**

Use the headballs in an upward direction from **bottom to top** along the jawline.

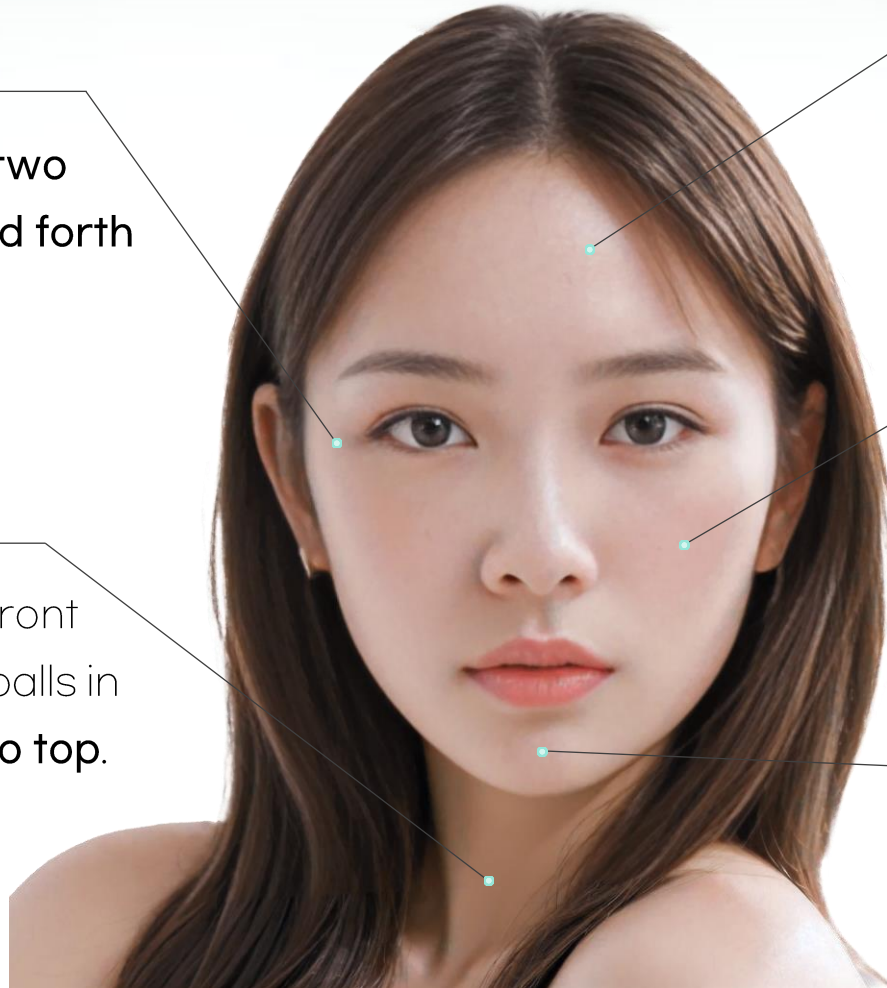

#HOW TO USE AURORA⁺



01

Please tidy up your skin after washing your face. Remove any remaining moisture on your skin.



02

Please apply your skincare products evenly.



03

Please activate the power button.



04

Please choose your desired level and use it on different areas for about 5 to 10 minutes.

Patents & Certificates

SAY SKIN



Title of Invention
“Atmospheric Pressure Plasma Device
Having Function of Providing Microcurrent”



Title of Invention
“Portable Beautifying Apparatus Using Plasma”



The background of the entire image is a flat-lay photograph of two white SAY SKIN Aurora skincare devices. The devices are sleek, handheld units with two spherical treatment heads at the top and control buttons on the side. They are surrounded by various skincare products, including jars of cream, tubes of gel, and a spray bottle, all arranged on a light-colored, textured surface. The overall aesthetic is clean and minimalist.

THANK YOU

Export Manager
: Chan Woo Jung

Email : chan@withnix.co.kr
Mobile : +82-10-3585-0464
Whatsapp : ChanWithnix
Wechat : ChanWithnix

Direct. 070-4338-3831
Withnix Co. Ltd /
<http://en.say-skin.com/>
<https://www.youtube.com/@sayskinbeautycare7628/>
© 2024 WITHNIX. All rights reserved