


How to keep vagina healthy with

we†trust®

Since 2003

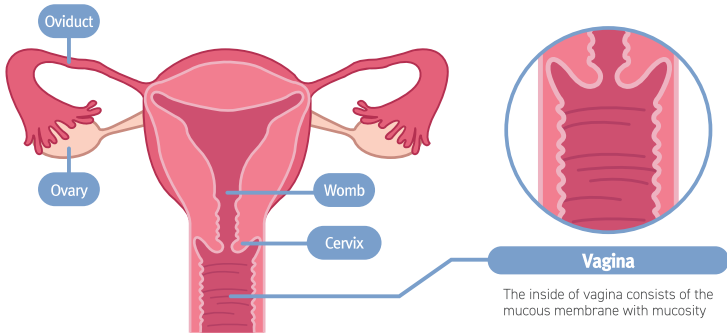

[VAGINA SELF-CLEANING ACTION]



Self-cleaning action of vagina has been proved by various researchers. According to various papers, glycogen (glucose) is the energy of the metabolism for Lactobacillus to produce Lactic acid. Lactic acid is one of the most crucial elements to keep the health of the vagina. Lactic acid supports the sub-acidic environment (pH 3.5-5.0) which promotes the growth of more Lactobacillus and the natural protective environment of the vagina. This is a simple cycle of how to keep a healthy vagina naturally itself.

[THE VAGINAL ECOSYSTEM]

Normal vaginal pH



Vaginal pH changes at various stages




[BY STAGES OF THE LIFE CYCLE]

	FEMALE HORMONE	LACTOBACILLI IN VAGINA	VAGINAL pH
Menarche & Puberty	⚙️ Unstable	▲ Appearing	💧 pH 3.8~4.5
Menstruation	↕️ Changing	⚙️ Unstable	⚙️ Unstable
Pregnancy	↕️ Changing	⚙️ Unstable	⚙️ Unstable
Menopause	▼ Slowly diminishing	▼ Slowly diminishing	🔴 pH 7.0

[BY STAGES OF THE MENSTRUATION]

- pH increase in period by menstruation blood (sub-alkaline).
- Increasing pH for surviving sperm in human physiology.
- Changing hormone by pregnancy, menopause and menstruation.

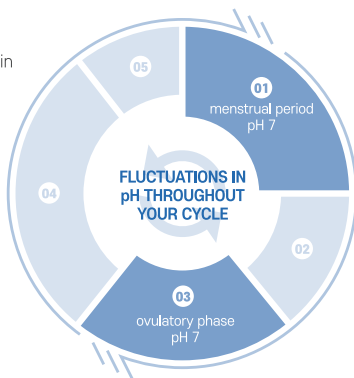

- 01** Blood pH ~7
Period 3~7 days
- 02** Normal Vaginal pH ~ 4.5
Follicular phase 11~12days
- 03** Cervical mucus pH ~7
Semen pH ~7
- 04** Normal Vaginal pH ~ 4.5
- 05** Luteal phase

FEEL HEALTHIER AND REFRESH WITH 'IN'CLEAR

By directly applying its cleansing gel of pH 4.2 ~ 4.8 into vagina.

As well known, healthy vagina is sub-acidic (an optimal balanced pH of 3.5-5.0) but, in case of douching or alkaline soap and body wash, vaginal flora can be disrupted and you are exposed to more risk on your vagina such as odor, itching, unusual discharge and vaginal inflammation. Therefore, the choice of right feminine care products requires a lot of careful attention. Now, wettrust confidently recommends you 'innovative & hygienic' disposable feminine cleanser of sub-acidic, water soluble and clear gel to use at any place & time.



Packing Unit :
1.7g*3pcs / 1.7g*10pcs

[DISPOSABLE FEMININE CLEANSER]

Medical device class 2 in Korea

Medical device class 2 in Japan

Help maintain vagina pH balanced with lactic acid

Individually wrapped portable size

Use without water at any place

Substantially effective at relieving odor, itching, irritation, abnormal discharge and preventing vaginal inflammation

CLEAN & MOISTURE

- ↳ Individually wrapped disposable applicator
- ↳ Water soluble transparent gel
- ↳ Containing Sodium hyaluronate & Lactic acid
- ↳ Cleansing gel (pH 4.2~4.8)



Sodium Hyaluronate



Lactic Acid




Sub-acid
(Cleansing gel)

Streamlined design for
smooth insertion into a vagina

Disposable applicator

Patented applicator of
ergonomic design

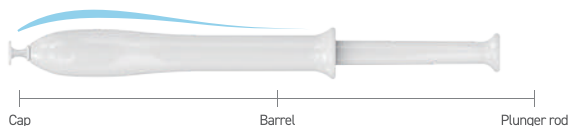

Containing Sodium hyaluronate
& Lactic acid



Cleansing gel 1.7g
Hypoallergenic transparent gel

Cleansing gel pH 4.2~4.8

INCLEAR
pH 4.2~4.8

GEL TYPE

Liquid

Semi-solid

Transparent

Opaque

TRANSLUCENCE



VAGINAL DRYNESS


(MENOPAUSE AND OTHER REASON)

When a woman stops her period specifically for lower estrogen, it is likely that there are itching, burning, irritation and dryness on her vagina since lower estrogen prevents vaginal wall from keeping healthy, thick, elastic and moisturized.

Now, vaginal dryness arises not only to women who are around and post menopause but also to a large number of women for all ages due to stress, life style, environmental and social factors.

[CAUSE OF VAGINAL DRYNESS]

During menopause, vaginal dryness is mainly caused by decreasing estrogen levels

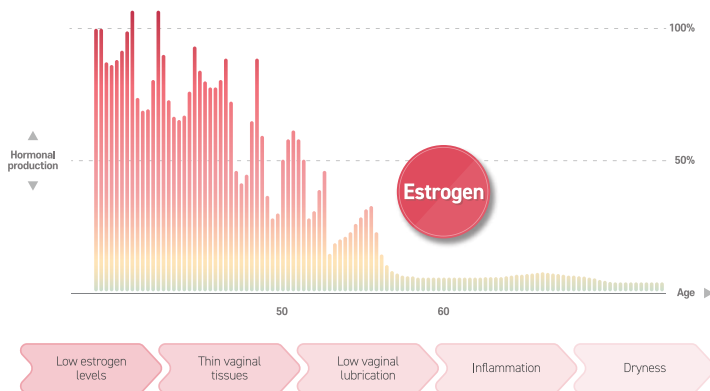

FEEL LIKE YOUR OWN NATURAL LUBRICATION !

wettrust VM is a pre-filled one-push syringe type vaginal moisturizer for short-term use. wettrust VM is designed for women who are experiencing vaginal dryness and it provides moisture beyond the menopause and immediate hydrating effects and comforts. wettrust VM provides relief from the symptoms of vaginal atrophy, vaginal dryness, and itching in women close to menopausal age.

Wettrust Vaginal Moisturizer

Packing Unit : 1.7g*3pcs / 1.7g*10pcs

Wettrust VM

Capacity : 50ml / 200ml



[VAGINAL MOISTURIZER]

Medical device class 2a in Europe

US FDA 510K cleared (Medical device class 2)

Supporting natural lubrication

Providing relief from the symptoms of menopausal and post-menopausal vaginal dryness

Sub-acid disposable vaginal moisturizer

Individually wrapped portable size

Use without water at any place

[APPLICATOR INSTRUCTIONS]



- 1) Wash hands and open the individual package using the easy-cut.
- 2) Take out approx. half of the applicator.
Hold the body of the applicator with your thumb and middle finger and place your index finger on the plunger lightly.
Then, hold the protective cap with fingers and remove it by bending at an angle of 45 degrees.

- 3) Take out the applicator from the individual package.
Get into a comfortable position and insert the applicator approx. 5-8cm. (from the tip)
- 4) Squeeze the gel into the vagina by gently pushing the plunger until it stops.
- 5) Pull out the applicator from the vagina and put the used applicator back into the individual package and throw away.

USE SAFELY, CONVENIENTLY!

Vaginal
Applicator

WETTRUST VA

- Easy-to-use device to effectively support insertion of various suppositories such as vaginal tablets.
- Relieving uncomfortable insertion of vaginal suppositories by hand.
- Streamlined ergonomic design for comfortable use and smooth insert into the vagina.
- Compatible with most of vaginal suppositories.



[The front of flower shape]

Soft and flexible material relieving a certain pressure women may feel whilst insertion into vagina.

Appropriately hold a vaginal suppository not to be out from the front.

* Vaginal suppositories are not included *



Round oval



Diamond



Bullet



Oval

Compatible with various vaginal suppositories under width 11.5mm

Packing Unit :
2pcs / 6pcs



[VAGINAL SUPPOSITORY APPLICATOR]

Medical device class 2 in Korea

Medical device class 1 in Japan

Medical device class 1 in Europe

Medical device class 1 in USA

Disposable suppository applicator

Easy to use for inserting various suppositories such as vaginal tablets

SIMPLE & SAFE

Wettrust Vaginal Applicator is ergonomic design for comfortable use and helps insert vaginal suppositories into vagina deep enough. This disposable applicator is individually wrapped and compatible with various vaginal suppositories.





Convenient & hygienic use



Simple & convenient
one push type



Individually wrapped

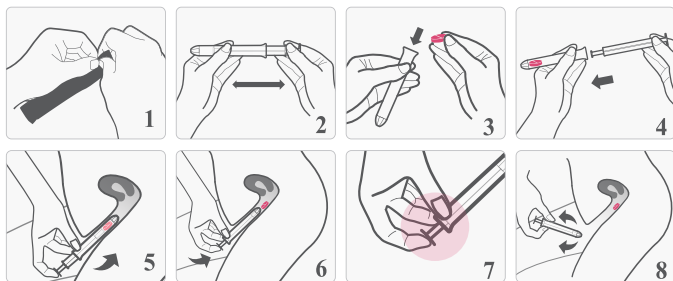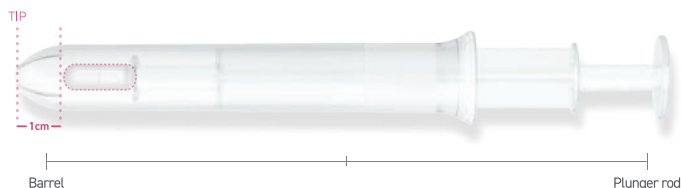

- 1) Wash your hands and remove the wrapper of wettrust Vaginal Suppository Applicator.
- 2) Make sure if the plunger moves back and forth smoothly.
- 3) Pull out the plunger and put a suppository (1cm from the tip) into the applicator.
- 4) Insert the plunger back into the applicator.
- 5) Take a comfortable position such as standing, squatting, bending forward at the waist, standing with one foot raised on a chair or stool.
- 6) Push the plunger to insert the suppository.
- 7) Make sure if the plunger does not move any longer so that the suppository settled in your vagina.
- 8) Remove the applicator completely and throw it away.



Website : www.wettrustglobal.com

E-mail : overseas@wettrustglobal.com

Tel : +82 · 31 · 719 · 9260 Fax : +82 · 31 · 719 · 8395

Address : 312 The Prau, 27 Jeongjail-ro, Bundang-gu,
Seongnam, Gyeonggi-do, Korea